

Micryo Plus® Cryogenic Labels



Dots and Strips

A new range of self-adhesive labels manufactured from matt white polyester, designed specifically for low temperature use. This dedicated range of labels can be used for most applications to temperatures as low as -196°C. They have a strong permanent adhesive (apply at +5°C or more). They can even be used in incubators.

The wide range of sizes allows you to label your whole experiment permanently from microcentrifuge tube to cryovial, microtitre plate to storage box.

Features:

Micryo labels can be used on polypropylene, glass, polycarbonate and most metals.

- No need to self-laminate
- Sleeker labels mean your tubes will still fit most storage systems (including cryocanes) without the label sticking or being torn off
- Easy to see, the bold range of 5 colours mean that you can identify your samples whilst still in the box
- Unique adhesive formulation for that extra staying power
- Guaranteed not to fall off under normal use (providing user instructions are followed)
- A cost effective solution for your low temperature labelling needs




For Tubes

PART NO.	COLOUR	TYPE	QUANTITY
----------	--------	------	----------




ASY/162	White	Dispenser Roll	Roll of 1000
ASY/163	White	Dispenser Roll	Roll of 5000
ASY/180	5 colours*	Dispenser Roll	Roll of 5000*
ASY/180R	Red	Dispenser Roll	Roll of 5000
ASY/180G	Green	Dispenser Roll	Roll of 5000
ASY/180P	Purple	Dispenser Roll	Roll of 5000
ASY/180Y	Yellow	Dispenser Roll	Roll of 5000
ASY/180B	Blue	Dispenser Roll	Roll of 5000

*(1000 of each Blue, Red, Green, Yellow, Purple)

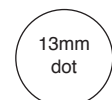

ASY/179	White	Sheets of 192 dots	Pack of 20 sheets
---------	-------	--------------------	-------------------

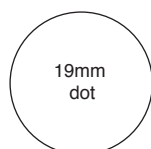

ASY/164	White	Dispenser Roll	Roll of 1000
ASY/165	White	Dispenser Roll	Roll of 5000

For Cryovials

Part no.	Colour	Type	Quantity
----------	--------	------	----------

20x6mm strip

ASY/161	White	Dispenser Roll	Roll of 1000 strips
---------	-------	----------------	---------------------

38x6mm strip

ASY/172	White	Sheets of 156 labels	Pack of 20 sheets
ASY/173	5 Colours**	Sheets of 156 labels	Pack of 25 sheets**

24x13mm strip

ASY/169	White	Sheets of 119 labels	Pack of 20 sheets
ASY/170	5 Colours**	Sheets of 119 labels	Pack of 25 sheets**

33x13mm strip

ASY/160	White	Dispenser Roll	Roll of 1000 strips
ASY/174	Blue	Dispenser Roll	Roll of 1000 strips
ASY/175	Red	Dispenser Roll	Roll of 1000 strips
ASY/176	Green	Dispenser Roll	Roll of 1000 strips
ASY/177	Yellow	Dispenser Roll	Roll of 1000 strips
ASY/178	Purple	Dispenser Roll	Roll of 1000 strips
ASY/167	White	Sheets of 85 labels	Pack of 20 sheets
ASY/168	5 Colours**	Sheets of 156 labels	Pack of 25 sheets**

** (5 sheets of each Blue, Red, Green, Yellow, Purple)

For Racks, Cryo and Freezer Boxes

67x25mm strip

ASY/171	White	Sheets of 30 labels	Pack of 20 sheets
---------	-------	---------------------	-------------------

Easy Remove Labels



Stick - Off Labels

Stick Off Labels

Stick-Off labels are ideal for applying onto laboratory glassware and plastics.

The benefits using this label is, you can either soak in warm water or wash off in one rinse cycle in a glasswasher. The label and adhesive bond will dissolve in warm water in less than a minute. Recommended temperature range for applied label use is -70°C to 60°C.

With six different sizes available in the range, the labels fit a variety of bottles, flasks beakers and other laboratory containers.

All labels are 100% biodegradable and can be written on with marker pen or pencil.



Printer Configuration

A comprehensive guide to sheet label page setup for Microsoft® Word® and Wordperfect® is available on request. Please call or email us for more details.

For Glassware and Plasticware

PART NO.	COLOUR	TYPE	QUANTITY
50x13mm			
ASY/5013	White	Dispenser Roll	Roll of 2500 strips
25x25mm			
50x25mm			
ASY/2525	White (25mm)	Dispenser Roll	Roll of 2500 strips
ASY/5025	White (50mm)	Dispenser Roll	Roll of 2500 strips
38x26mm			
ASY/3826	White	Dispenser Roll	Roll of 2500 strips
50x38mm			
ASY/5038	White	Dispenser Roll	Roll of 2500 strips
76x32mm			
ASY/7632	5 Colours**	Dispenser Roll	Roll of 2500 strips

Custom Labels

Our labels are manufactured by us and can be produced in customised sizes, shapes and materials according to your requirements. Please note that pricing for this customised label service is quantity dependant.

29 Hall Barn Road Industrial Estate
 Isleham
 Ely, Cambs
 CB7 5RJ
 Email: sales@asynt.com
 Tel: 01638 781709