
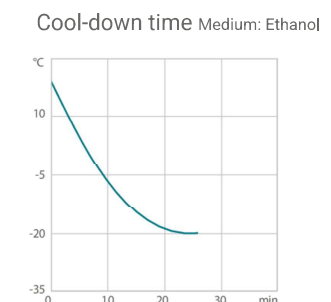
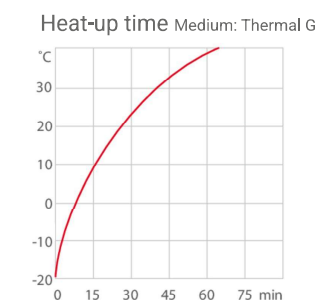


Product data sheet

VALEGRO 500 Recirculating cooler

Modern, user-friendly chillers for a wide range of cooling tasks

The VALEGRO 500 recirculating cooler with a working temperature range -20 ... +40 °C offers, in relation to its very compact design, a strong cooling capacity of 0.5 kW at 20 °C. The integrated pump generates a pump pressure up to 0.85 bar. It is also infinitely adjustable and can be ideally adapted to the requirements of the customer application. The modern units offer precise temperature control for a wide range of applications in laboratories, industry and research and fit seamlessly into any working environment thanks to their compact design. The refrigeration circuit works with an environmentally-friendly, natural refrigerant. When designing the VALEGRO recirculating coolers, the focus was on maximum user-friendliness. Many helpful features such as a timer function, a removable venting grid that can be removed without tools or an illuminated and simple to clean level indicator make daily work easier. The inclined OLED touch display is easy to read and ergonomic to operate. The intuitive menu navigation provides simple and quick access to all functions. Warnings and messages in easy-to-understand full-text language ensure safe handling of the units at all times. In addition to RS232 and USB-C interfaces, the VALEGRO recirculating coolers are also available with Ethernet (optional). All interfaces are located directly under the display and are easy to use from the front.

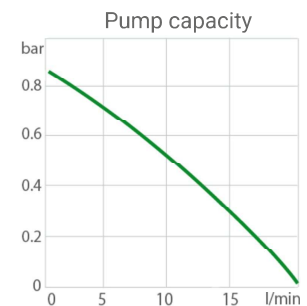

Product Features

- ✓ Precise temperature stability of ± 0.3 °C
- ✓ Adjustable pump for regulating pressure and flow rate
- ✓ Removable venting grid without tools
- ✓ Integrated castors for good handling on and under the lab bench
- ✓ Simple timer function
- ✓ Autostart function
- ✓ Data recording
- ✓ Intuitive menu navigation in full text German/English



Available options

Interface

Ethernet
VALEGRO Ethernet activation



Performance values

230V/50Hz (UK Plug Type BS1363A)

Viscosity max. cSt	50
Pump capacity flow pressure l/min	19
Pump capacity pressure bar	0.85
Power consumption A	3
Voltage Tolerance %	230V ±5

Order number including voltage version (230V/50Hz) 9610050.30003286

Cooling capacity (Ethanol)

°C	20	0	-10	-20
kW ¹	0.5	0.37	0.25	0.1

Note about natural refrigerants:

Temperature control units using natural refrigerants are often subject to regulatory requirements regarding the installation site, operation, transport or disposal of the units. If you have any questions, we will be happy to advise you.

Refrigerant stage 1

Refrigerant R290

Filling weight g 46

Global Warming Potential for R290 3

Carbon dioxide equivalent t 0.000138

¹ Performance specifications measured in accordance with DIN 12876. Cooling capacities up to 20 °C measured with ethanol; over 20 °C with thermal oil unless otherwise specified. Performance specifications apply at an ambient temperature of 20 °C. Performance values may differ with other bath fluids.

Technical data

Available voltage versions

Order No. 9610050

Available voltage versions:

9610050.30003285	100V 50/60Hz (Nema N5-15 Plug) (R290)
9610050.30003284	115V 60Hz (Nema N5-15 Plug) (R290)
9610050.30003283	200-230V 50/60Hz (Schuko Plug - CEE 7/4 Plug Type F) (R290)
9610050.30003286	230V 50Hz (UK Plug Type BS1363A) (R290)
9610050.30003287	230V 50/60Hz (CH Plug Type SEV 1011) (R290)
9610050.30003288	200-230V 50/60Hz (CN Plug) (R290)

Cooling

Cooling of compressor 1-stage Air

Bath

Bath tank Stainless steel

Electronics

Interfaces	Ethernet optional, RS232, USB
Temperature control	PID1
Temperature display	OLED
Temperature setting	Touchscreen
Electronic Timer h:min	99 ... 59

Dimensions and volumes

Weight kg	26
Barbed fittings inner diameter mm	8/12 mm
Total dimensions cm (W x L x H)	23.5 x 45 x 45
Body dimensions cm (W x L x H)	23 x 40 x 44.5
Filling volume l	1.5 ... 3
Pump connections	M16x1 male

Temperature values

Return flow temperature max. °C	+85
Working temperature range °C	-20 ... +40
Temperature stability °C	±0.3
Ambient temperature °C	+5 ... +40
Temperature display resolution °C	0.1

Other

Sound pressure level dbA	55
Classification	Classification I (NFL)
Standards	CE, RoHs
Alarm notification	Acoustic, Optical
IP Code	IP 21
Pump function	Pressure Pump
Pump type	Centrifugal Pump
User Interface Language	English, German



Included in delivery

2 Barbed fittings for tubing 8 and 12 mm ID. (Pump connections
M16x1 male)

Refer to www.julabo.com for more information regarding the entire JULABO product portfolio. Technical changes without prior notification.
Images may deviate from the original. | Datasheet No.EN9610050/250729